



## AMERICAN INTERNATIONAL UNIVERSITY – BANGLADESH

Department of Natural Science (Chemistry)

Faculty of Science & Technology

Programs: B.Sc. Eng'g (EEE/CSE/IPE)

**CHEM 1101: CHEMISTRY**

# Chemistry Lab Report

**Semester:** Spring

**Session:** 2022-2023

**NO EXPERIMENT, NO REPORT**

**Experiment No:**  7

**Name of the Experiment:** Determination of ferrous ion ( $Fe^{2+}$ ) in a supplied solution of iron salt by standard potassium dichromate ( $K_2Cr_2O_7$ ) solution.

**Date of Performance:** 29-03-23

**Date of Submission:** 05-04-23

**Course-Teacher:** Dr. Mohammad Tariqul Islam

### **Instructions:**

1. A lab report consists of three parts: a cover page, body of the report and a data and results sheet (lab-sheet).
2. This is the cover page of a report and students will collect and preserve the lab-sheet of a particular experiment to be performed.
3. Body of the report includes-(1) Objective of the Experiment, (2) Theory, (3) Name of the Chemicals, (4) Name of the Apparatus, (5) Percentage of Error (if necessary) and (6) Discussion (I. Precautions taken, II. Possible errors).
4. Use A4-size off-set paper, write on one side of the paper by hand keeping suitable margin.
5. Staple the lab-sheet at the end of the report and cover page on the top.
6. Submit the report in time to avoid deduction of marks.
7. Students working in a group will write and submit the report individually.
8. Copying of the report from others is strictly prohibited.

**Name of the Student:** TRIDIB SARKAR  
**ID No:** 22-46444-1      **Section:** F      **Group:** 03

**FOR FACULTY USE ONLY**




**Faculty comments:** ..... , **Signature:** .....  
**Date:** .....

Objective: To know the amount of iron ( $\text{Fe}^{2+}$ ) in a supplied solution of iron salt by standard  $\text{K}_2\text{Cr}_2\text{O}_7$  solution.


Theory:

(i) Method: A redox titration is based on the neutralization reaction between the analyte and titrant. This one used a redox indicator to determine the endpoint. Frequently either the reactants or the titrants have colour intense enough that an additional indicator is not needed.

(ii) Reaction: Mohr's salt is ferrous ammonium sulphate to determine the amount of iron, a known value of solution titrated with  $\text{K}_2\text{Cr}_2\text{O}_7$  with diluted  $\text{H}_2\text{SO}_4$ .  $\text{K}_2\text{Cr}_2\text{O}_7$  oxidises  $\text{FeSO}_4$  present in Mohr's salt into  $[\text{Fe}_2(\text{SO}_4)_3]$ .

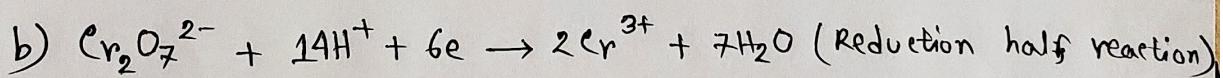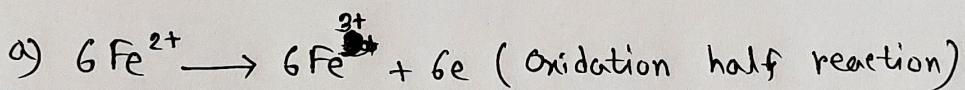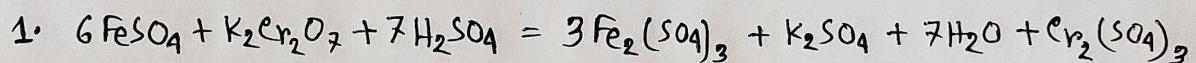

(iii) Indicator: Diphenyl amine is the organic compound with the formula  $(C_6H_5)_2NH$ . It is a colourless solid, but samples are often yellow due to oxidized impurities. It shows a very clear colour change from green to violet when end point of the titration is reached.

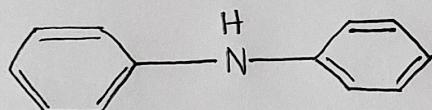

Diphenyl amine,  $(C_6H_5)_2NH$

Required Chemicals:

Name of the chemicals	Chemical formula
1. Iron salt solution	$\text{FeSO}_4 \cdot (\text{NH}_4)_2\text{SO}_4 \cdot 6\text{H}_2\text{O}$
2. 5% Sulfuric acid	$\text{H}_2\text{SO}_4$
3. Cone. Phosphoric acid	$\text{H}_3\text{PO}_4$
4. Standard Potassium Dichromate solution	$\text{K}_2\text{Cr}_2\text{O}_7$
5. Diphenyl Amine indicator	$\text{Ph}_2\text{NH}$

Apparatus:

- |                                |                          |
|--------------------------------|--------------------------|
| (i) Burette (50 ml)            | (v) Watch glass .        |
| (ii) Pipette (10 ml)           | (vi) Pipette filler .    |
| (iii) Conical flask (250 ml)   | (vii) Dropper .          |
| (iv) Volumetric flask (100 ml) | (viii) Stand and clamp . |

Name: TRIDIB SARKAR ID No: 22-46444-1 Section (Group): F (o3)

AMERICAN  
INTERNATIONAL  
UNIVERSITY –  
BANGLADESH (AIUB)



## Experiment 7


CHEM 1101: CHEMISTRY (EEE/CoE/CSE/IPE)

### EXPERIMENT NO. 7: DETERMINATION OF FERROUS ION ( $\text{Fe}^{2+}$ ) IN A SUPPLIED SOLUTION OF IRON SALT BY STANDARD POTASSIUM DICHROMATE ( $\text{K}_2\text{Cr}_2\text{O}_7$ ) SOLUTION.

**OBJECTIVE:** To know the amount of iron ( $\text{Fe}^{2+}$ ) in a supplied solution of iron salt by standard  $\text{K}_2\text{Cr}_2\text{O}_7$  solution.

#### **THEORY:**

- (i) Method: Redox titration  
(ii) Reaction:

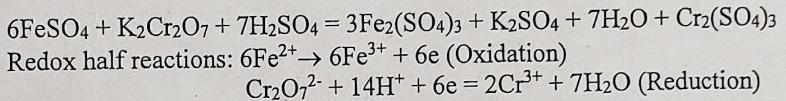

- (iii) Indicator: Diphenyl amine,  $(\text{C}_6\text{H}_5)_2\text{NH}$

#### **APPARATUS:**

Burette (50mL), pipette (10mL), conical flask (250mL), volumetric flask (100mL), watch glass, pipette filler, dropper, Stand and clamp etc.

#### **REQUIRED CHEMICALS:**

1. Iron salt solution,
2. 5% Sulfuric acid,
3. Conc. Phosphoric acid,
4. Standard  $\text{K}_2\text{Cr}_2\text{O}_7$  solution,
5. Diphenyl Amine indicator

#### **PREPARATION OF APPROX. 0.1N POTASSIUM DICHROMATE SOLUTION.**

Transfer approx. 0.49 gram of pure  $\text{K}_2\text{Cr}_2\text{O}_7$  into a 100 mL measuring flask and then dissolve it with distilled water up to the mark.

$$\text{Strength of } \text{K}_2\text{Cr}_2\text{O}_7 \text{ solution (S)} = \frac{\frac{0.51}{\text{Weight taken (in gm)}} \times 0.1}{0.49} (\text{N})$$

$$= 0.10408 \text{ N}$$

(Expt. 7 contd.)

**PROCEDURE:** Take 10 mL of the supplied iron salt (Mohr's salt) solution in a conical flask. Add 50 mL 5% sulfuric acid and 5 mL of conc. phosphoric acid. Then add 4-5 drops of diphenyl amine indicator and titrate slowly against the standard potassium dichromate solution drop wise maintaining an interval of few seconds between each drop until the addition of one drop causes the formation of intense purple or violet blue coloration which remains permanent and is unaffected by further addition of dichromate solution. Repeat the experiment at least thrice. Calculate the amount of iron per 500 mL of iron salt solution.

### EXPERIMENTAL DATA:

**Table:** Determination of the amount of iron in Mohr's salt solution using standard  $K_2Cr_2O_7$  solution.

No. of reading	Vol. of Mohr's salt solution (in mL)	Vol. of $K_2Cr_2O_7$ (burette reading) (in mL)			Mean (in $\mu\text{L}$ ) (V)
		Initial	Final	Difference	
1	10	0	9.70	9.70	
2	10	4.70	9.50	4.80	
3	10	9.50	14.50	5.00	
4	10				4.83

**CALCULATIONS:** 1 mL 1N  $K_2Cr_2O_7$   $\equiv$  0.05584 gm of  $Fe^{2+}$

$$\begin{aligned} \text{Amount of iron in 10 mL of iron salt solution} &= 0.05584 \times V \times S \text{ gm} \\ &= 0.05584 \times 4.83 \times 0.10408 \text{ gm} \\ &= 0.02807 \text{ gm} \end{aligned}$$

$$\begin{aligned} \text{Amount of iron in 500 mL of iron salt solution} &= 0.05584 \times V \times S \times 50 \text{ gm} \\ &= 0.05584 \times 4.83 \times 0.10408 \times 50 \text{ gm} \\ &= 1.40356 \text{ gm} \end{aligned}$$

$$\text{Observe value of } Fe^{2+} \text{ (in 500mL solution)} = 1.40356 \text{ gm}$$

$$\text{Known value of } Fe^{2+} \text{ (in 500mL solution)} = \frac{55.84 \times 10}{392.14} = 1.42398 \text{ gm}$$

### RESULTS:

Amount of  $Fe^{2+}$  ions in 500 mL of supplied Mohr's salts solution is  
1.40 gm

### PERCENTAGE OF ERROR:

$$\frac{\text{Known value} - \text{Observed value}}{\text{Known value}} \times 100 = \frac{1.42398 - 1.40356}{1.42398} \times 100 = 1.43\%$$

### Students should know

- Why it is necessary to use both the *sulfuric acid* as well as *phosphoric acid* in the reaction?
- Atomic weight, molecular weight of  $K_2Cr_2O_7$  and  $KMnO_4$ .
- Could you use  $KMnO_4$  instead of  $K_2Cr_2O_7$ ?
- Why the solution shows light bottle green colour after addition of  $K_2Cr_2O_7$ .

Discussion :

(i) Precautions:

- a) The apparatus was washed first two times with normal water then one time with distilled water.
- b) Care was taken when adding  $H_2SO_4$  in the solution.
- c) Care was taken while adding  $H_3PO_4$  solution.

(ii) Possible errors:

- a) Error might be occurred when taking the burette reading.
- b) Error might be occurred when observing the value of  $Fe^{2+}$  in 500 ml solution.
- c) Might have added more than 4-6 drops of  $(C_6H_5)_2NH$ .
- d) Concentration and calculation error might be occurred.